
NAIPO
FIND WELL-BEING

User Manual

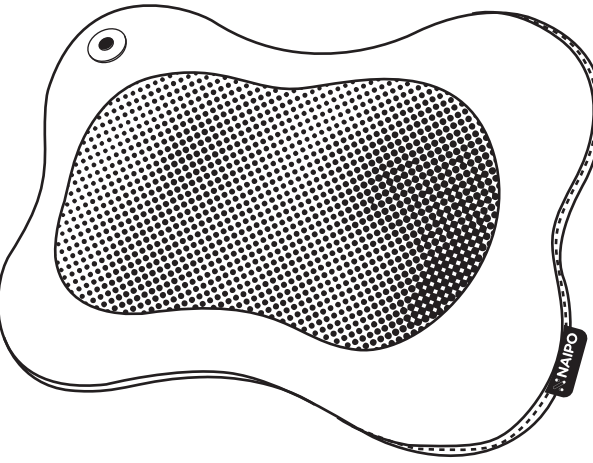Shiatsu Pillow Massager with Heat for Back & Neck
Model: MGP-129M

Shenzhen Pinking E-Commerce Co., Ltd.
Web: www.naipocare.com
E-mail: support@naipocare.com
Add: Huanan City Electronic Trading Building P07
Room 120, Longgang District, Shenzhen, China.

CE FC RoHS
Made in China

Contents

| | |
|----------|-------|
| English | 01-05 |
| Deutsch | 06-10 |
| Français | 11-15 |
| Español | 16-20 |
| Italiano | 21-25 |

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Important Safety Instruction

When using the product, all instructions and basic precautions should always be followed, including the following:

DANGER

To reduce the risk of electric shock: Always unplug this product from the electrical outlet immediately after using and before cleaning.

WARNING

To reduce the risk of burns, fire, electric shock, or injury to persons:

- The product should never be left unattended when plugged in. Unplug from outlet when it is not used and before putting on or taking off parts.
- Do not operate under a blanket or a pillow. Excessive heat can occur and cause fire, electric shock, or injury to persons.
- Do not carry this product by its supply cord or use the cord as a handle.
- Never operate the product with the air openings blocked. Keep the air openings free of lint, hair, and the like.
- Do not drop or insert any object into any opening.
- Do not use outdoors.
- Do not operate where aerosol (spray) products are being used or where oxygen is administered.
- To disconnect, switch off all controls and remove the plug from outlet.
- Close supervision is necessary when this appliance is used by, or near children, invalids, or disabled persons to ensure that they do not misuse the product.

- 01 -

- Keep the cord away from heated surface, otherwise it may cause serious burns. Do not use over insensitive skin areas or in the presence of poor circulation.
- The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.
- The appliance must only be supplied at safety extra low voltage core sponding to the marking on the appliance.
- The appliance is to be used with the power supply unit provided with the appliance.
- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or persons lacking experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Special attention should be given to components most susceptible to wear and it's advised to replace defective components immediately and/or keep the equipment out of use until repair.
- Use this appliance only for its intended use as described in this manual.
- Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or if it has dropped into water. Return the appliance to a service center for examination and repair.

- 02 -

The device is not to be used for commercial use or medical purposes. If you have health concerns, consult your doctor before using the product.

Only use indoors!

Do not use the device in damp rooms (e.g. when bathing or showering)

• Do not wash!
• Do not dry clean!
• After using this product for 20 minutes, rest for 30 minutes, then use again.


SAVE THESE INSTRUCTIONS

- 03 -

Massage Areas

Product Diagram

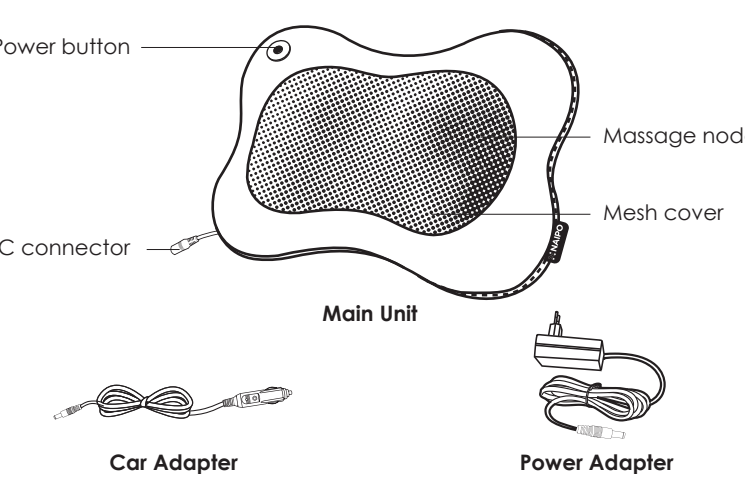Product Operation Guide

- Connect the product to the power supply.
- Press the power button to turn on the product, with kneading massage and heating function being activated.
- Press the power button again to turn off the heating function.
- Long press the power button for 3 seconds to stop all functions.

After using, disconnect the product from the power supply.

- 04 -

Product Maintenance

- Before cleaning, please make sure that the product is off and that it's disconnected from the power supply. Please clean using dry cloth.
- Keep the product clean and far away from humid environment.

Attention

Please make sure that the product is not cleaned using benzene, thinner, or other aggressive detergents, otherwise the product color might fade.

Contact Us

If you need any assistance or support, please contact us. We will do our best to provide you with a viable solution.
Web: www.naipocare.com
Email: support@naipocare.com

Technology Parameters

Name: Massage Pillow
Model: MGP-129M
Rated Voltage: DC 12V ± 2.5A
Rated Power: 30W

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